



# Single Use Hemorrhoidal Circular Stapler

**Tissue Anti-Retreat Mechanism Design**  
32mm Lower than PPH03 (4mm)

- Deeper and Larger Staple Housing
- Easier to Retain More Cutting Tissue



## Major Features

- Accommodate to II~IV internal hemorrhoids
- Accommodate STARR procedure to treat Outlet Obstructed Constipation/Syndrome
- Protect anal cushion and retain precise faeces control
- Simple and minimum invasive to use
- Faster procedure, less pain, quicker recovery

**Golden Stapler Surgical Co., Ltd.**

CE ISO CFDA

0434 9001/13485:2003

*Advancing Surgical Technology  
to Benefit More People*

# Single Use Hemorrhoidal Circular Stapler

## Staple Housing is 6mm Longer than PPH03

Staple housing size is larger than PPH03 over 30%  
Staple housing is able to contain much more cutting tissue



## Tissue Anti-Retreat Mechanism Design

32mm Lower than PPH03 (4mm)  
Much easier to retain more cutting tissue into staple housing



## Plastic Cutting Ring

Optimized tactical & audible responses during firing



## Less Distance between Firing and Grabbing Handle

Facilitate female surgeons' easier operation  
Silicon handle avoids slippage and maintains a great hand touch feeling



## Transparent Purse-String Suture Anoscope

Optimum surgery view and position chosen during Purse-String suturing



## Larger Suture Threader Ports

Significantly facilitate dragging cutting tissue into staple housing



## Blunt Design of the Staple Housing Edge

Easier to drag cutting tissue into staple housing  
Avoid tissue damage during dragging



## Butterfly Shape Anal Dilator

Accommodate to obese, strong hip muscle or local anesthetic patients




## Magnified Green Indication Window with a Gap Setting Scale

Controllable staple formation height between 0.8-2.0 mm  
Accommodate to PPH and STARR procedures  
Effectively achieve suturing and homeostasis purposes



## Non-detachable Central Rod Design

One-piece central rod design to avoid anvil fall off  
Solid and non-detachable anvil offers premium staple formation basis

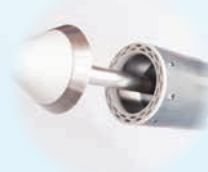

## Size 32 Contains 4 More Staples than PPH03

Better off anastomosis and homeostasis



Code	Open Leg Length (MM)	Staple Qty	Closed Staple Height (MM)	Outer Dia. (MM)	Circular Knife Dia. (MM)	Staple Rows
PPH32	4.5	32	0.8-2.0	32.5	24	2
PPH34	4.5	34	0.8-2.0	34.5	26	2

**Golden Stapler Surgical Co., Ltd.**

Add: Building 7A, Jiangsu Wujin Sci - Tech Pioneer Park, 256 Mingxin Zhong Road,  
Hutang Town, Wujing District, Changzhou, Jiangsu Province, China  
Tel: +86 519 85172266 Fax: +86 519 85172299  
Web: <http://www.goldenstapler.com> P.C: 213164