



# Single Patient Use Linear Cutter and Subassembly

## Patented Design of Curved Cutting Knife

- Sharper cutting & easier firing



### Major Features

- Transection of GI tract
- Side-to-side anastomosis of GI tract
- Tailoring & suturing of gastric tissue
- Tailoring & suturing of pulmonary tissue
- Suturing of GI tract incision

Golden Stapler Surgical Co., Ltd.

CE ISO SFDA

0434 9001/13485:2003

Advancing Surgical Technology  
to Benefit More People

### Intermediate Position

One-handed manipulation  
Repositioning on targeted tissue  
Precise cutting & suturing of targeted tissue



### Tissue Retaining Pin

Avoid tissue slippage out of instrument distal end  
Ensure effective cutting length




### Parallel Closure

Even compression on tissue  
Uniform staple formation height  
Effective suturing & hemostasis



### Protruding Cam Mechanism

Adopt lever principle  
Establish parallel closure

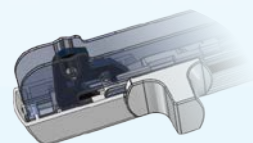

### Curved Cutting Knife Design

Shaper cutting & easier firing



### Blue, Yellow & Green Cartridge Reloads

Appropriate to various targeted tissue thickness  
Satisfied suturing & hemostasis



### Suturing Rows longer than Cutting Slot

Effective suturing of distal cutting tissue  
Avoid bleeding of distal cut tissue



### Scale Measurement

Clear indication of side-to-side anastomosis length



Code	Colour	Staple Rows	Staple Qty	Closed Staple Height (MM)	Suturing Length (MM)	Applicable Tissue Thickness	Code of Subassembly Cartridge Reload
LCA55W	White ●	4	56	1.0	55	Blood vessel, thin tissue	LCAR55W
LCA55B	Blue ●	4	56	1.5	55	Normal thick tissue	LCAR55B
LCA55Y	Yellow ●	4	56	1.8	55	Medium thick tissue	LCAR55Y
LCA55G	Green ●	4	56	2.0	55	Thick tissue	LCAR55G
LCA75B	Blue ●	4	76	1.5	75	Normal thick tissue	LCAR75B
LCA75Y	Yellow ●	4	76	1.8	75	Medium thick tissue	LCAR75Y
LCA75G	Green ●	4	76	2.0	75	Thick tissue	LCAR75G

Golden Stapler Surgical Co., Ltd.

Add: Building 7A, Hutang Kechuangyuna 256 Mingxin Zhong Road,  
Hutang Town Wujin, Changzhou, Jiangsu, China

Tel: +86 519 85172266

Fax: +86 519 85172299

Web: <http://www.goldenstapler.com>

P.C: 213164